



P-808®

Hair removal- Faster, painless
and better

CONTENTS



- 1 P-808 perspective
- 2 Working Principle
- 3 Operator interface
- 4 Technical Advantages
- 5 System Applications
- 6 Specifications
- 6 System advantages

P-808 perspective



Super cooling

Sapphire Cooling 0-5°C

Professional

Micro channel technology

High

power Germany bar

Huge

12*12mm /12*23mm



HANDPIECE

- Modern Design
- Stable output
- Short time irradiation
- Permanent effect
- Easy to control
- Excellent clinical effect
- High security

Working Principle

- P-808 System, use special laser with long Pulse-Width 808nm, can penetrate to hair follicle .
- Using selective light absorption theory laser can be preferentially absorbed by heating the hair shaft and hair follicle, moreover to destroy the hair follicle and oxygen organization around hair follicle.
- When laser outputs, system with special cooling technology, cool the skin and protect skin from being hurt and reach a very safe and comfortable treatment.



Normal Hair Follicle



Instant Hair Follicle after Treatment

Golden yellow: the light of confidence, healthy skin

Turquoise blue: pure, clear, energetic.

Black: mystic, majestic, full of attraction.

Pink: gentle, romantic, filled with illusion.

Gray: independent, undertone, symbol of mature




Operator interface


P-808 Hair Laser Removal System




P-808

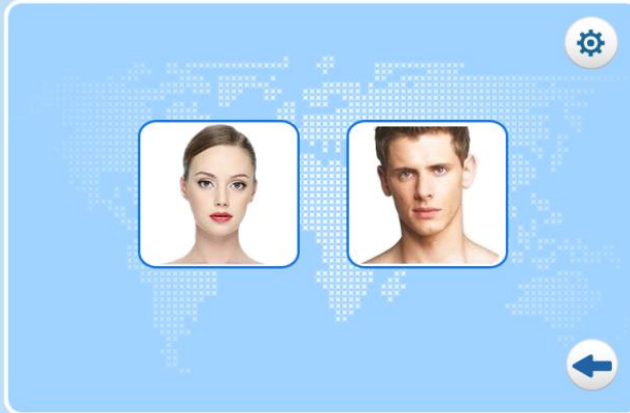

P-808

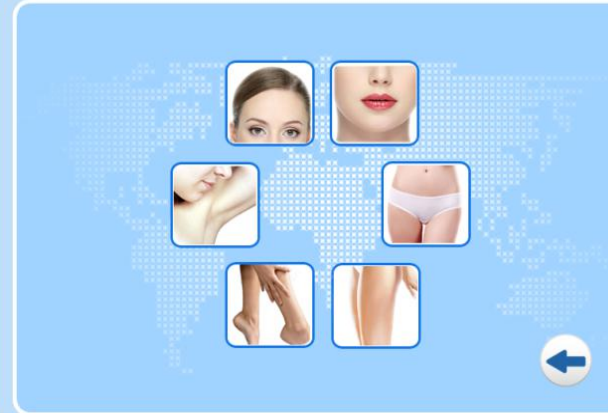

P-808

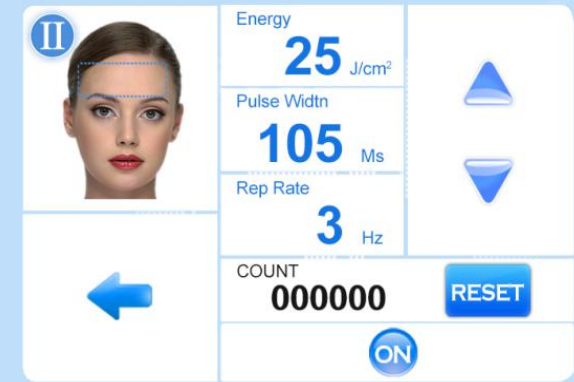

1. Touch screen is the main operation and display part of the hair removal system.
2. Adopting the blue and white color, the interface is concise, readable and considerate.
3. According to the ergonomics, the operator interface is designed separately for both male and female.
4. Through the touch screen, you can achieve the functions include shutdown, laser parameter settings and ready and use the laser.
5. You could only do the parameter setting, mode switch, refrigerator opening and shutdown during standby mode.

Specifications

Wavelength	808nm
Spot size	12*12mm, 12*23mm
Repetition rate	0.5-10HZ
Pulse width	10-1200ms
Laser output	600W
Energy density	1-120J/cm ²
Sapphire contact cooling	-5~0°C
Operate interface	12' True color touch screen
Cooling system	air & water cooling
Power supply	85-260V full range,50/60HZ
Laser Bar	10 bars
Life time of laser bar	10 million shots
Dimension	
Weight	35Kgs

Technical advantages

- High Energy Output: Make sure good effective for any coloring matter -----120J
- Long Pulse Width: Make sure heating hair follicle, lasting hair removal-----1200ms
- Big Spot Size: Fast and Effective Results -----12*12/12*23
- Strong Cooling: Instance cooling skin, safe and comfortable ----- -5°C
- Fast Speed: slide treatment to make it safety and convenience-----10HZ

System Applications

- Treat all kind of hair color from black hair to white hair;
- Treat all skin types from white to dark skin;
- No pain and shorter treatment sessions ;
- Effective and safe treatment for permanent hair removal;



System advantages

1. Long laser pulse width

Laser pulse duration is 10ms to 1200ms adjustable, can do treatment effectively special for different hair and lasting!

2. Super big true color touch screen make the operation more easy

12 inch multi color touch screen, you can choose man or woman , skin color and treatment areas,very easy to operate !



System advantages

3. No pigmentation, used for any skin type.

P-808 hair removal system uses penetrate technology and skin absorbs fewer lasers and no pigmentation.


4. Pure Sapphire contact cooling system can make sure that treatment more safe and no hurt to skin.

Pure Sapphire contact cooling system of P-808 hair removal system can make sure that doctor use the maximum energy and can not hurt skin. This kind of cooling system can make sure that no hurt to skin when patients are given more pressure. As reported, the pressure can increase 30%-40% light.



System advantages

5. The big spot size can speedup the treatment speed.
P-808 use 12*12mm big spot size. It makes sure that when do large area hair removal and prompt treatment speed and efficiency. Because of big spot size, it can escape from repeated treatment for hair area.
6. The power module design is suitable for import and export business. It is convenient to maintenance.
7. Just like the picture show right.

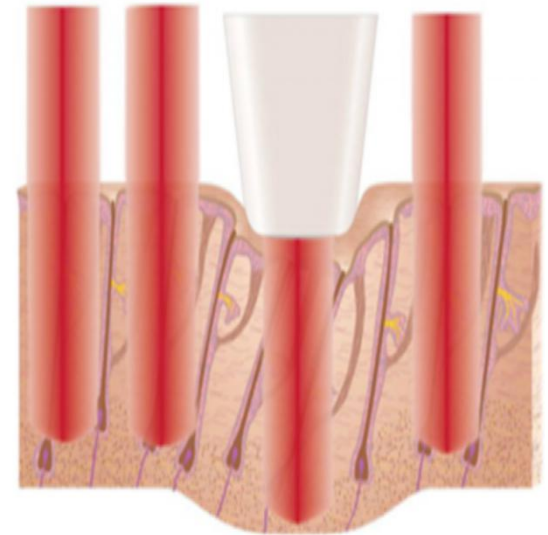

Before and After Treatment





Established in 2005, Beijing sano Beauty S & T Co., Ltd takes the leading in applying advanced laser and Intense light technology in medical and beauty industry.

Now we developed Mona IPL series, E-light, Co2 fractional laser, 808nm diode laser, ultrasonic cavitations, fractional RF and hair re-growth laser.

With efficient R & D team, mass production ability and reliable quality control system, Beijing Sanhe Beauty supplies the OEM & ODM service and help our clients build up their own reputation and more competitive in the market.

All the process of the producing is under ISO13485 Quality system and match with CE certification request.

Beijing Sanhe Beauty will work hard constantly for supplying best product and after sale service.

Real Picture of Beijing Sanhe Beauty S & T Co.,Ltd

THANKS

Welcome to consult more!

www.sanobeauty.com

Address: 2nd floor, Building No.7, Maohuagongchang, Nancaizhen, Shunyi District,

BeijingTel: 86-15810588766, 0086-53006037

Email: gxl@sanhebeauty.com

